

VARIO®化学隔膜泵 ME 16C NT VARIO select

[网站链接](#)

产品说明

对于连续无油泵抽腐蚀性气体和蒸汽，ME 16C NT VARIO select泵是一个很好的解决方案。该单级泵性能优异，抽速大，极限真空度达70mbar，可用来蒸发大量的低沸点溶剂。采用耐化学腐蚀氟塑料的泵头盖，以及经过验证的聚四氟乙烯三明治夹层结构隔膜，保证了我们的化学隔膜泵具有无与伦比的使用寿命和维护间隔。集成式的VACUU·SELECT真空控制器搭配了简单易用、基于应用程序的用户界面，内嵌了所有常见的实验室真空应用。VACUU·SELECT真空控制器可满足您所有要求。对于简易实验过程可手动修改控制参数，也可进行全自动蒸馏，或者通过简单的“拖放”操作编辑创建自己的应用程序。在溶剂蒸发过程中，VACUU·SELECT真空控制器能够持续探测溶剂沸点，并自动调节真空过程，维持过程控制。通过变频控制，泵只按需要的速度运行。这延长了维修间隔并且泵运行安静。

性能特点

- 配置图形用户界面及预定义应用的VACUU·SELECT真空控制器能简化实验室工作
- 自动沸点探测及真空度调节能够缩短反应时间
- 适用于大量低沸点溶剂的蒸发
- 超长的膜片寿命，维护周期长且维护成本低

订购信息

产品名称	额定电源电压 工作频率	电源插头类型	订货号
ME 16C NT VARIO select	200-230 V 50-60 Hz	CEE	25741750

货品配置

泵体安装完毕，连接即可使用，附带说明书


配件

- VACUU·SELECT**用户管理** VACUU·SELECT User Management (20700900)
- DN 8-10mm (20636588)
- **小法兰的中心支架** KF DN 20/25 (20660001)
- **小法兰外齐式环形中心支架** KF DN 20/25 (20660196)
- **PVC真空管，内置不锈钢螺线，带小法兰** KF DN 25 (20686021)
- **带软管喷嘴接头的小法兰** KF DN 25 / DN 10 mm (20662807)
- **橡胶真空管** DN 10 mm (20686002)
- **密封圈套装** ME 8C NT (20696867)


技术指标

50 Hz 最大抽速 (mbar)	19.3 m³/h
60 Hz 最大抽速 (cfm)	11.4 cfm
极限真空	70 mbar / 53 torr
极限真空, 使用气镇	100 mbar / 75 torr
泵头数	8
泵级数	1
真空控制器	VACUU-SELECT
的允许环境温度范围 (运行)	10 - 40 °C
的允许环境温度范围 (存放)	-10 - 60 °C
最大出气口压力 (绝对压力)	1.1 bar
进气口连接	KF DN 25 小法兰
出气口连接	DN 15 mm 软管喷嘴接头
额定电机功率	0.53 kW
电机转速范围	30 - 2400 min ⁻¹
防护等级	IP 40
尺寸 (长x宽x高)	533 mm x 260 mm x 450 mm
重量	28.9 kg
噪音水平 50 Hz/1500 Upm/62% (VARIO)/1500 Upm (VARIO-SP)/12500 Upm (VACUU-PURE®)	47 dBA
ATEX规格	II 3/- G Ex h IIC T3 Gc X Internal Atm. only
NRTL 标准	否
通讯	Ethernet (RJ45): Modbus TCP / USB (Type A): Data transfer, software update, RS-232 (additional adapter required)


50 Hz下的抽速曲线

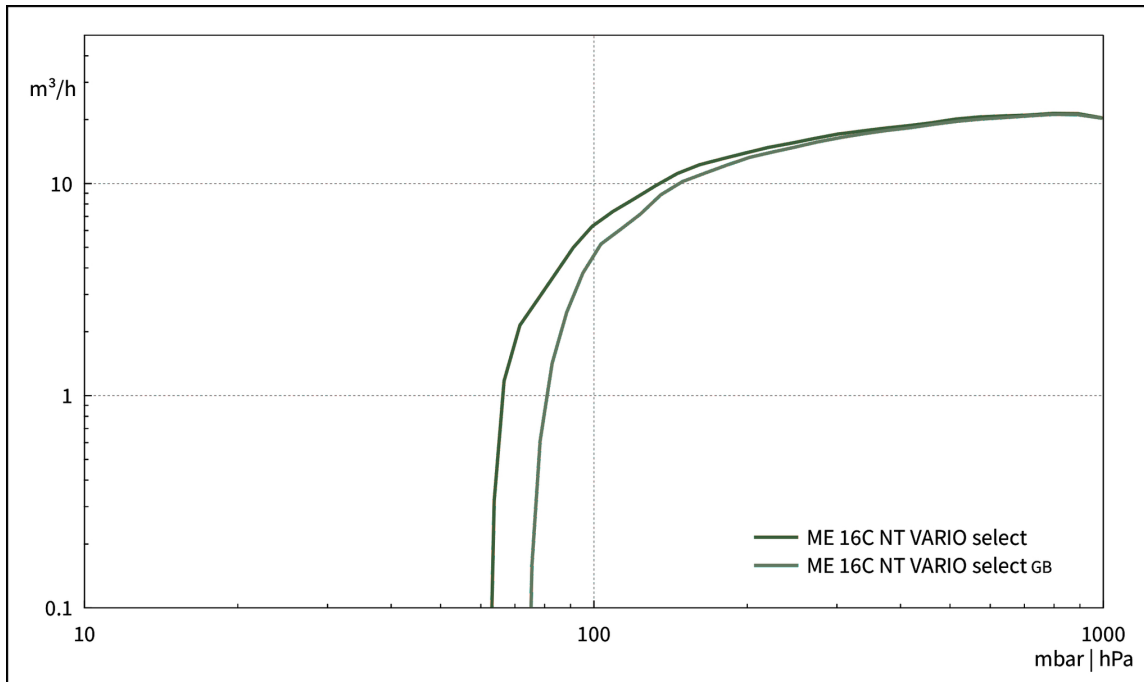


50 Hz下的抽气曲线（100升容积）

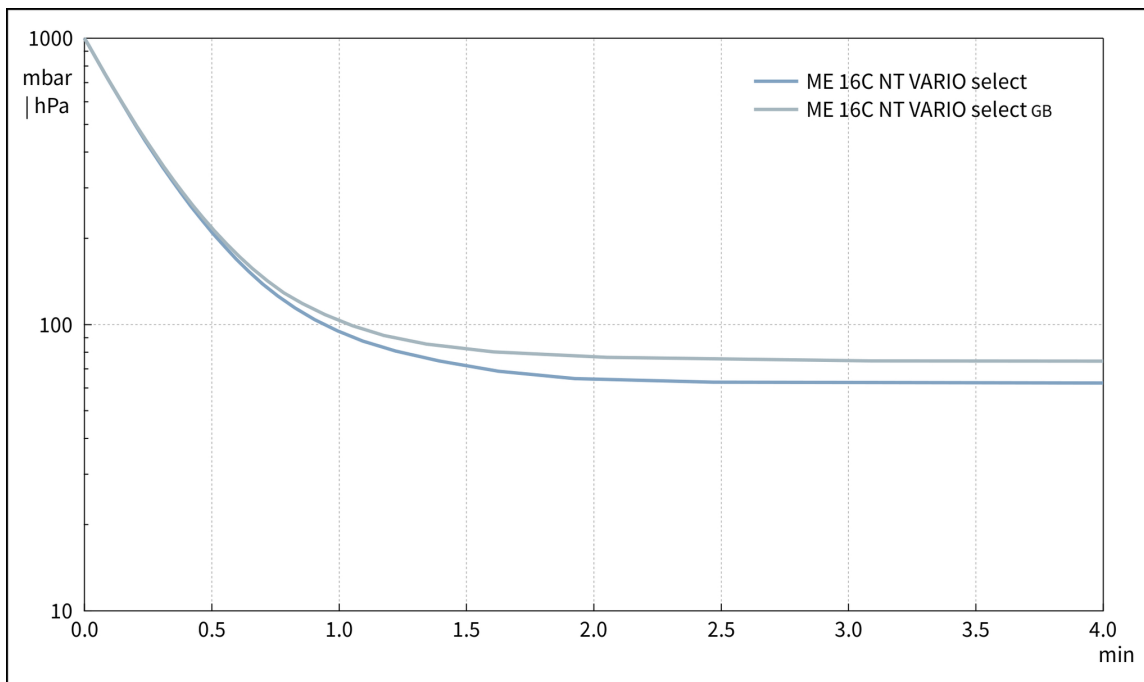


60 Hz下的抽速曲线

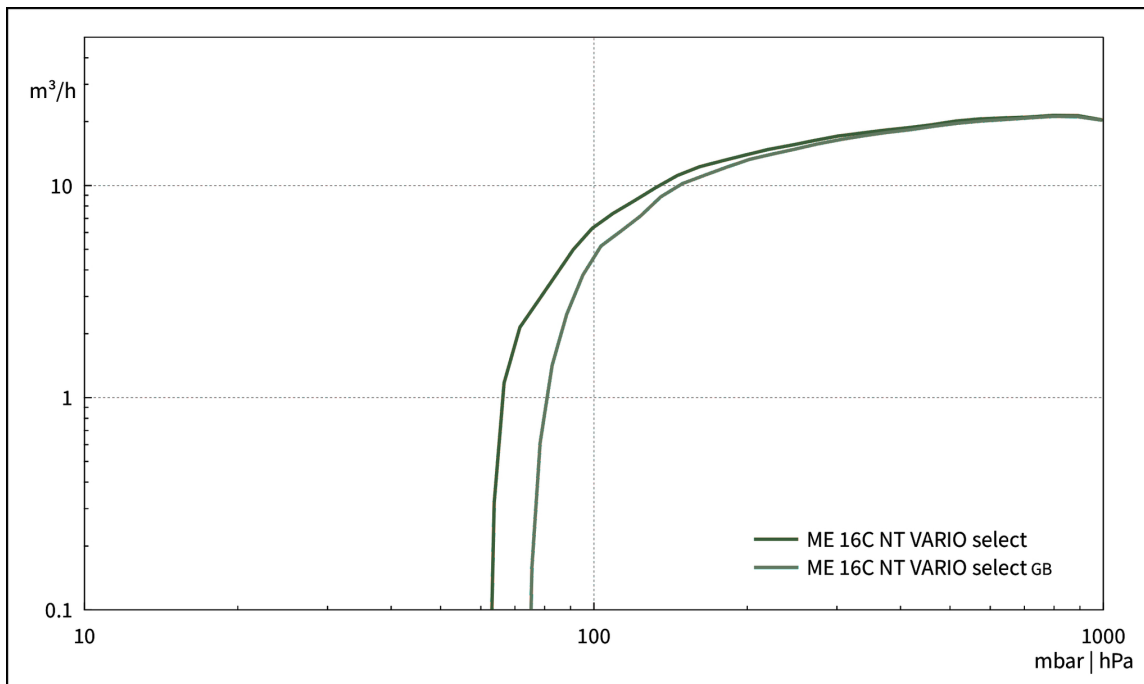

60 Hz下的抽气曲线（100升容积）

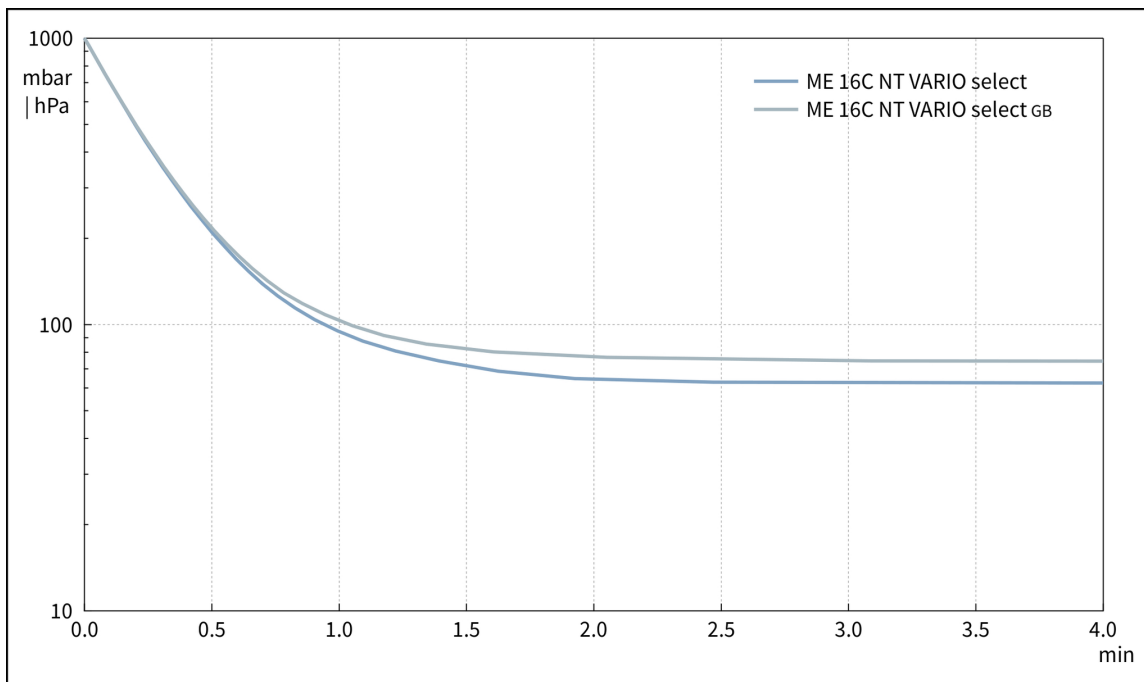

尺寸规格表

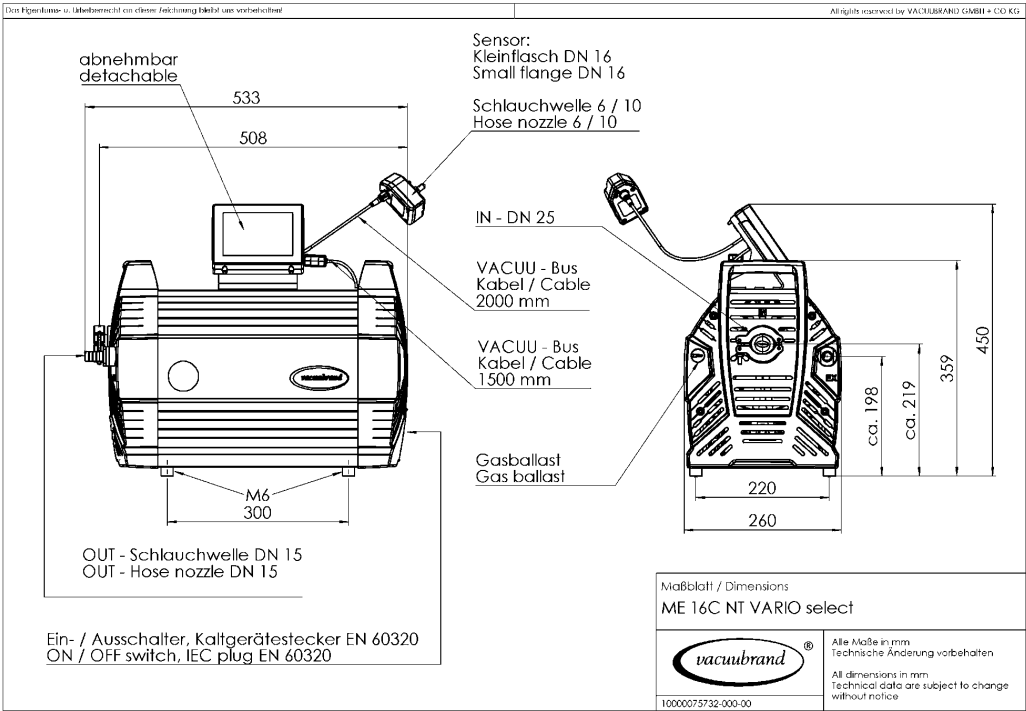

VACUUBRAND GMBH + CO KG
P.O. Box 1664
97877 Wertheim
T +49 9342 808 5550
F +49 9342 808 5555

info@vacuubrand.com
www.vacuubrand.com

VACUUBRAND®

VACUUBRAND is part of the Brand group, a manufacturer of high-quality and innovative laboratory equipment and vacuum technology with a global team of approx. 1,000 employees.

With highly motivated and qualified employees, we manufacture the world's most comprehensive product family for vacuum generation, measurement and control for the rough and fine vacuum range in the laboratory.

All with one goal:

Enabling our customers to run reliable, safe and efficient processes. This is why we commit ourselves to continuous innovation, first-class performance, and outstanding quality.

Engineering excellence is central to our ability to deliver products and solutions renowned for their quality, durability, and reliability. Innovative technology, high quality materials and extensive testing deliver unmatched performance - engineered and manufactured in Germany.

Together with you, our experts will find the best solution for your application and provide you with first-class service along the way.

© VACUUBRAND GMBH + CO KG

BRAND (Shanghai) Trading Co., Ltd.
Shanghai, China

T +86 21 6422 2318
info@brand.com.cn
china.brand.com.cn

BRAND Scientific Equipment Pvt. Ltd.
Mumbai, India

T +91 22 42957790
customersupport@brand.co.in
www.brand.co.in

BRANDTECH® Scientific, Inc.
Essex, CT. United States of America

T +1 860 767 2562
info@brandtech.com
www.brandtech.com